Türkçe / English

KK-06 03.2021 / Rev:00

KULLANIM KILAVUZU

Instructions for Use

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Read instructions before use!

Warning!

CE marking according to Directive 93/42 EEC, referring to the manufacturer's notified body

Non-Sterile

Store in a dry place

Protect from sunlight

Lot/Batch Number: Located on both the product and the package.

Catalog Number: Located on both the product and the package.

Production Date: Located on the package.

INTENDED USE AND INDICATIONS

Aspiration and Irrigation cannulas are used to make the operation area suitable by aspirating fluid accumulations, blood and tissues during surgery, to remove an obstruction or to irrigating solutions and medicines into the operation area.

CONTRAINDICATIONS

- Uterine cannulas are contraindicated for pregnant patients and may cause miscarriage.

The Aspiration and Irrigation Instruments are contraindicated for all other applications except their intended use! The conditions that can risk the success of the operation such as injuries or accompanying diseases should be evaluated by the user. Risks should be taken into account for anatomical structures in the planned area.

GENERAL SAFETY WARNINGS

- Read, the instructions for use before each use.
- Keep this instructions for use throughout the product use.

- Products are **NON-STERILE!** It must be cleaned and sterilized before use according to the instructions below.

- Used only by people the requisite training, knowledge and experience.
- Do not use out of its intended use.
- Before each use, instruments should be checked / inspected. (See **Inspection & Check**)
- Do not use out of sight or in situations where vision is not available.
- In cases with infectious disease risk, the legislation determined by the relevant authority institution must be followed.
- Before each use, check the intercompatibility of instruments that are used with another device.

WARNINGS REGARDING CLEANING, DISINFECTION AND STERILIZATION

- For the cleaning, disinfection and sterilization of products, national and international standards, directives and clinical hygiene instructions must be followed.
- Use cleaning and disinfectant agents according to the manufacturer's instructions.
- Apply cleaning processes (manual or mechanical) to newly purchased instruments before the first sterilization process.
- For patients with a diagnosis or suspicion of Creutzfeldt-Jakob disease (CJD) or variants, national regulations regarding disinfection and sterilization of products must be followed.
- Handle with validated methods and calibrated autoclaves and machines at all times.
- Use Personal Protective Equipment (PPE) when working with contaminated instruments.PPE: apron, mask, glasses, gloves, shoe covers etc.
- Prefer mechanical cleaning cleaning instead of manual cleaning as it gives more reliable results.
- Handle with care. Do not drop or apply excessive forces to the instruments during reprocessing.

TECHNICAL DESCRIPTION

Structure / Shape

1- Shaft
2- Luer Hub
3- Finger Cut-off*
4- Stylet**


*Finger cut-off provides control of the suction pressure. Available on some models.

**The stylet helps to clean the residues that cause blockage/obstruction in the cannula. Available on some models.

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- Aspiration and irrigation cannulas work by connecting to aspiration/irrigation devices with a suitable hose connection from the Luer hubs.

Suction and Irrigation Cannula **(Standard)**

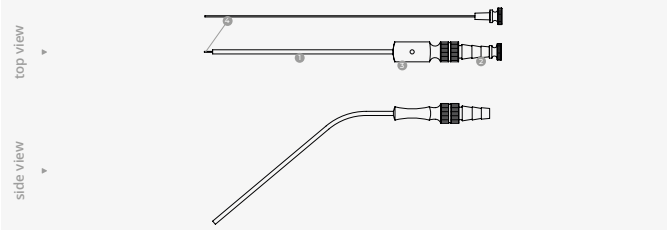Cleaning and Maintenance Products / Tools

warning **The use of metal brushes or clothes cause corrosion!**

warning **The use of chlorine solutions (saline solutions, bleaching agent, ringer solution) cause pitting corrosion or stress corrosion cracking!**

- Use cleaning and disinfection agents appropriate for aluminum, high-grade steel.
- Use cleaning agents with approved efficacy (VAH or FDA approval or CE mark)
- Cleaners with a disinfectant effect (manual cleaning) must be in accordance with EN 14885 or equivalent national guidelines.
- Use soft, lint-free fabrics, paper cloth, plastic, soft bristle brushes and water spray guns
- Do not use aldehyde based cleaning products due to its protein fixing properties.
- Cleaning solutions should be changed on daily basis. If the solution gets contamined visibly, it can be changed earlier.
- If powder products are used, the powder should be completely dissolved in water before use. Undissolved particles may cause changes in instrument surfaces and clog narrow tube instruments.
- Compressed air drying is more effective than drying with fabric or cloth.
- Ultrasonic cleaning is recommended to contribute to both manual and mechanical cleaning especiallyfor lumened and articulated instruments
- Suitable detergent should be selected for ultrasonic cleaning. (Non-foaming)
- Separate stainless-steel instruments from other metal instruments during ultrasonic cleaning to avoid electrolysis.
- Be sure that products are free of cleaning agent at the final rinse of ultrasonic cleaning step.

Water qualifications

Softened water*: Water hardening substances (calcium and magnesium cations) are reduced.

Demineralised water*: Mineral substances largely removed by reverse osmosis or cation & anion exchangers or electrode ionization or distillation.

For steam sterilization and the last rinsing of cleaning processing, limit values for feed water quality is specified in EN 285 and ISO 17665.

⚠ Non-conforming water properties cause:

- Scaling, lime deposits,
- Brown-red deposits, secondary rust,
- White-grey, colored appearance, thin scaling,
- Pitting,
- Spotting and scaling.

* See www.ak-i.org, Publications, Red Brochure for water properties (softened water, demineralized water)

CLEANING, DISINFECTION AND STERILIZATION INSTRUCTIONS

Point of Use

- Remove excess biologic soil from the instruments with a disposable wipe.
- Rinse non-visible surfaces with demineralised water, with a disposable syringe.
- Place dry products in sealed waste containers.

Containment and Transportation

⚠ warning — Wet disposal conditions increase the risk corrosion and microbial growth. The residues of cleaning agents interfere with the following processes.

- Transport the instruments to decontamination are in closed containers at dry state.
- Keeping the products at a dirty state more than 6 hours may cause corrosion and fail cleaning processes. Used instruments should be processed within 6 hours at dry disposal.
- Prevent slipping of the instruments by using special carrying tools during transportation.

Preparation for Cleaning / Precleaning

⚠ warning — Do not exceed 25 °C for pre-cleaning temperature!

- Rinse the instruments to remove coarse contaminations.
- Disassemble the detachable instruments to preparing for cleaning.
- Ultrasonic washing is recommended before disinfection to increase the cleaning efficiency of heavily soiled or difficult to clean instruments such as lumened/tubular, micro and articulated instruments for mechanical support.
- If a pre-cleaning step is necessary before automatic cleaning/ disinfection, follow the 1st and 2nd process steps in the "**Manual Cleaning/Disinfection**" section of this document. Then send the instruments to the automatic washing-disinfection process by dry transfer.

Manual Cleaning / Disinfection

⚠ warning — Do not exceed the maximum permitted cleaning temperature of 45 °C. Otherwise protein denaturation may ocur.

⚠ warning — Use a non-protein fixing cleaning solution does not contain aldehyde, phenol and QUAT with a disinfection effect according to the manufacturer's instructions.

⚠ warning — All scrubbing should be performed below the surface of the solution. Do not apply manuel cleaning to the cutting instruments or sharp pointed instruments to avoid physical demage to the operator !

	Process Step	Water Quality	Heat	Time
1	Cleaning/disinfecting[1]	Softened water*	Room Temp.	15 min
2	First Rinsing	Softened water*	Room Temp.	1 min
3	Disinfecting[1]	Softened water*	Room Temp.	15 min
4	Final Rinsing	Demineralized Water*	Room Temp.	1 min
5	Drying	-	Room Temp.	**

1. Pour a 25 ml dose cleaning/disinfectant solution in 5 liters of water (%0,5 dilution). Completely immerse instruments in the solution and gently shake them to remove trapped bubbles. Actuate instruments with hinges or moving parts to ensure contact of the solution with all surfaces. Flush the solution into non-visible surfaces and holes of lumened and tubular instruments such as cannulas with a disposable syringe (20 ml). While soaking scrub surfaces using a suitable cleaning brush until all visible soil has been removed. Actuate moveable mechanisms. Particular attention should be given to hinged joints, box locks, instrument teeth, rough surfaces and areas with moving components or springs. Lumens, blind holes and cannulas should be cleaned using a snug fitting round suitable brush. Insert the brush into the lumen, blind hole or cannula with a twisting motion while pushing in and out multiple times.

2. Remove the instruments from the solution and rinse with softened water for a minimum of 1 minute. Actuate all moveable and hinged parts while rinsing. Thoroughly and aggressively flush lumens, holes, cannulas and other difficult to access areas. Allow excess water to drip of from the instrument prior to disinfectant solution.

3. Pour a 125 ml disinfectant solution dose in 5 liters of water (%2,5 dilution). Immerse fully in the disinfectant solution. Move the joints at least 5 times beneath the solution. Flush the solution into lumens with a disposable syringe (20 ml).

4. Rinse with running demineralized water for at least 10 seconds close and open the joints at least 5 times, flush lumens with a disposable syringe. We recommend using fully desalinated water (demineralized, drinking water quality at least) to avoid water spots.

5. Dry instruments with a clean, absorbent non-shedding lint free cloth or compressed air to remove moisture from lumens, holes, cannulas and difficult to access areas.

[1] Anios - Aniosyme DD1 Solution was used to verify cleaning.

- Check for cleanliness visually at the end the cleaning process especially the hinges, lumens, ratchets, teeth and etc. Repeat the cleaning process if necessary.
- Apply maintenance products (see **Maintenance & Care**)
- Control visually and functioning of the instruments according to the "**Inspection & Check**" part of this document.

Mechanical cleaning and thermal disinfecting

warning

Use low alkaline cleaners (8.0-11.0 pH). Do not use high alkaline value (>11 pH) cleaning agents.

warning

The washer disinfecfer must have either FDA approval or CE mark (validated according to EN ISO 15883)

- Place instruments in a washer/disinfector. Follow the washer / disinfector manufacturer's instructions for loading the instruments for maximum cleaning exposure; e.g. open all instruments, place concave instruments on their side or upside down, use baskets and trays designed for washers, place heavier instruments on the bottom of trays and baskets. If the washer/disinfector is equipped with special racks (e.g. for cannulated instruments) use them according to the manufacturer's instructions.
- Arrange the position of the instruments in the washer/disinfector to prevent mechanical damage through contact.
- Avoid overloading to make sure the surface of the instrument effectively gets in contact with cleaning solution.
- To prevent corrosion, unload the washer/disinfector immediately when the cleaning program ended.
- The following minimum wash cycle parameters are validated and recommended

	Process Step	Water Quality	Heat	Time
1	Pre-wash	Softened water*	<25 °C	3 min
2	Cleaning	Demineralized Water*	40 – 55 °C	10 min
3	First rinse / Neutralization	Demineralized Water*	>10 °C	2 min
4	Intermediate rinse	Demineralized Water*	>10 °C	1 min
5	Thermal disinfection / Final rinse	Demineralized Water*	90 °C	5 min
6	Drying	-	-	**

- Check for cleanliness visually at the end the cleaning process especially the hinges, lumens, ratchets, teeth and etc. Repeat the cleaning process if necessary.
- Apply maintenance products (see **Maintenance & Care**)
- Control visually and functioning of the instruments according to the **"Inspection & Check"** part of this document.

***** See "Water qualifications" section of this document for water properties.

****** Dry time is dependent upon the load size placed into the washer/disinfector.

Maintenance & Care

warning

Be sure to apply care product (lubricant) to metal sliding surfaces, moving parts, joints/box locks prior to **"Inspection & Check"**.

warning

Raw materials used in care products must comply with pharmacopeias such as USP/NF, Ph. Eur, and national and international cosmetics regulations. **Do not use** silicone oil based products!

- Ensure that the instruments are cleaned effectively before maintenance. Visually and physically inspect the instruments for blood and protein soils. It is recommended to verify the cleaning with a chemical protein test, especially on surfaces that are difficult to clean or invisible such as ratchets, hinges, teeth, lumens.

Care products must be:

- Paraffin/white oil based,
- Biocompatible,
- Suitable for steam sterilization and vapor-permeable.

- Allow the products to cool down to room temperature.
- Apply the care product directly to joints, hinges and friction surfaces by using a squirt oiler, an oil pen, or through targeted application of drops of oil.
- Actuate moving mechanisms a few times to distribute the care product evenly.
- Wipe the excessive care product with a lint-free cloth.

warning

Failure to apply sufficient care products to moving parts causes friction related damage and corrosion!

Inspection & Check

Verify the following:

- Check the laser marking of the product. It must be easily readable visually or mechanically.
- There must not be any crack, fracture, wearing, burr, discoloration, spotting, rust or corrosion on the surface of the product.
- There must not be any loose, bent, broken, cracked, worn, or fractured components.
- Lumened/tubular instruments must be checked against blockages with a suitable cleaning brush. If a blockage is detected, reprocessing must be repeated.
- Matching parts must be compatible with each other for the assembled products.

warning

If any of the above controls are not confirmed, do not use the product, contact with Bahadır Medical Devices (manufacturer).
Damaged products must only be maintained by the manufacturer.

Packaging prior to Sterilization

- Pack the products according to the determined sterile packaging. *****
- Ensure that all cutting edges and fine tips are protected.

*** The Bahadır Sterilization Containers were used in the validated sterilization method.**

Sterilization Method

Moist heat/steam sterilization is the preferred and recommended method for Bahadır Reusable Surgical Instruments.

- The steam sterilizer must comply with the EN 285 standard and be verified according to the EN ISO 17665 standard.
- Ensure that the maximum permitted load according to the manufacturer's instructions is not exceeded.
- For steam sterilization, limit values for feed water quality as specified in EN 285 and ISO 17665.
- Ethylene oxide or other sterilization methods are not recommended.
- The use of chemical indicators causes stains on the instrument surface.

Sterilization Method	Heat	Time
Prevacuum	134 °C[1,3]	5 min[2,3]

[1] Validated exposure time required to achieve a 10-6 sterility assurance level (SAL)

[2] Validated exposure temperature required to achieve a 10-6 sterility assurance level (SAL).

[3] Local or national specifications should be followed where steam sterilization requirements are stricter or more conservative than those listed in this table.

Storage

- Storage and distribution of the product should be maintained by the designated staff according to the rules of the instutition which is used.
- Newly purchased products should be stored in a dust and moist-free environment.
- Surgical instruments should not be stored near to chemicals.
- Sterile, packaged instruments should be stored in a designated, limited access area that is well ventilated and provides protection from dust, moisture, insects, vermin, and temperature/humidity extremes.

Life Span

- As long as the suction and irrigation cannulas are regularly maintained, used according to their intended use, storage and maintenance conditions are followed, they can serve for 7 years depending on the average usage intensity.
- Products which can not pass from functional and visual inspection can not be used.

Disposal

- The disposal of a product that comes into contact with the human body and its packaging should be done in accordance with the national laws of the country where the product is used..

Technical Service

- Guarantee term: 2 years.
- The user is responsible for any adverse situation that may arise from the use of products whose warranty expires.

Manufacturer Contact Information

BAHADIR TIBBİ ALET CİHAZ VE İNŞAAT MAKİNA SANAYİ VE TİCARET ANONİM ŞİRKETİ
Şabanoğlu OSB Mah. Vali M.Erdoğan Cebeci Bul. No:29/1
Tekkeköy - Samsun / Turkey

Tel: +90 362 266 70 00
Fax: +90 362 266 70 05
info@bahadir.com.tr / www.bahadir.com.tr

Notified Body Contact Information

KIWA BELGELENDİRME HİZMETLERİ AŞ
(İTOSB) İstanbul Tuzla Organize Sanayi Bölgesi
Tepeören Mevkii 34957
Tuzla - İstanbul / Turkey

Tel: +90 (216) 593 25 75
posta@kiwa.com